


# User Manual

# Content


Product Structure .....	02
Product Instruction .....	03-04
Product Specification .....	04
How To Use Device .....	05-06
Applicable Skin Color and Hair Color .....	06
Applicable Parts and Recommended Use Time .....	07
Common Problems and Solutions .....	07-09
QA .....	09-10

# Product Structure

## ● Product form

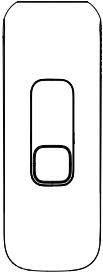


Main Machine

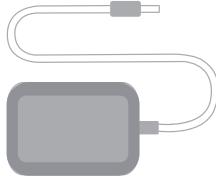

Adapter





User Manual

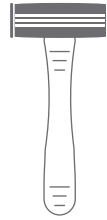


Razor

## ● Product structure

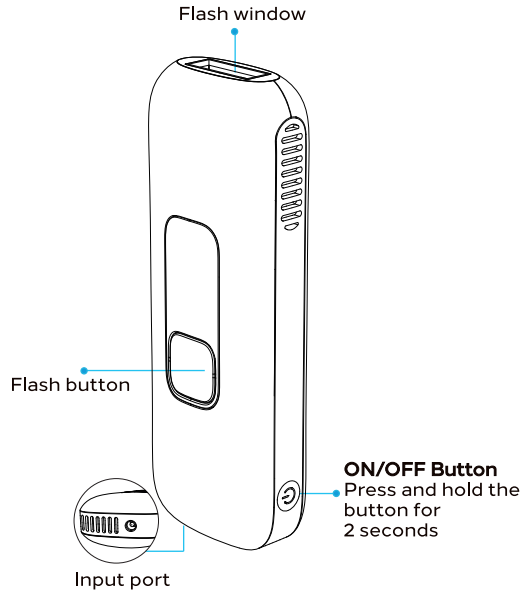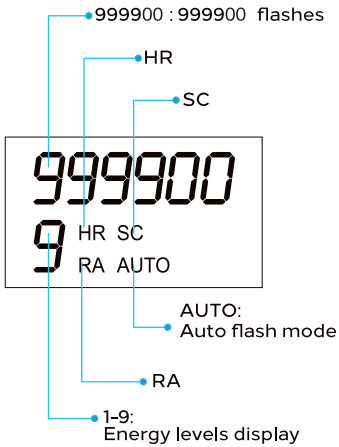

## Product Instruction

Under the irradiation of pulsed light of a specific wavelength and pulse width, the melanin in the skin absorbs the energy and transmits the transient heat to the hair follicle. When the hair follicle reaches a certain temperature, it enters a dormant phase, causing the hair to stop growing in stages and eventually fall out. Finally, by using the device for about 8 weeks or more, the desired results can be achieved.

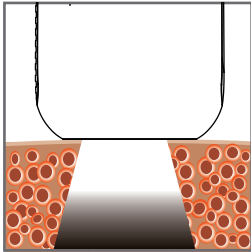

### ● Product introduction

Multi-functional hair removal device is an innovative beauty product that effectively and safely removes unwanted body hair. It uses intelligent pulsed light irradiation technology and combines with optical technology for beauty care, which is one of the most professional and practical products in this field

### ● Please read the manual carefully before use


The ultimate goal of Desianina hair removal device is to ensure the comfort and safety of the client. Please read all operating instructions, safety precautions and warnings carefully before use.

### ● How IPL works

With Intense Pulsed Light technology, gentle pulses of light are applied to the skin and absorbed by the hair root. The lighter the skin and the darker the hair, the better the pulses of light are absorbed.


The pulses of light stimulate the hair follicle to go into a resting phase. As a consequence, the hair sheds naturally and hair regrowth is prevented.

● **Use 8 weeks as a complete treatment.**

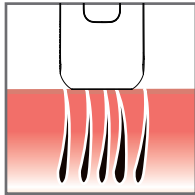


**Start (1-4 weeks)**

Treat twice a week for the first 1-4 weeks.

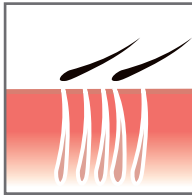

**In between (5-8 weeks)**

Weekly treatment during weeks 5 to 8

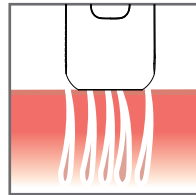

**Maintenance (After 8 weeks)**

Continue IPL to your needs (use every 1-2 months).


## Product Specification

Product name:	Hand-held IPL device
Adapter input and output:	100-240V~50/60Hz 24V/ 1.5A
Light source:	IPL pulse light
Energy levels adjustment	9 levels adjustable
Working humidity:	the relate humidity 30%-60%
Light output wave band:	600-1200nm
Light output area:	3 cm <sup>2</sup>
working temperature:	41-82°F / 5-28°C
Product size:	179*64*39mm
Maintain guarantee limit period:	1 year

## How To Use




1. Shave the hair before each treatment to get better result.

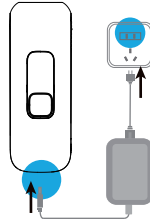

2. Connect power supply.




3. Wear Goggles to protect your eyes. ( Note: Do not look directly at the flash.)

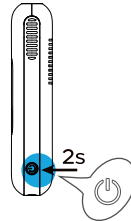


4. Long press the "power" button for 2 seconds to turn on the device, and the energy indicator will light up. (Machine startup defaults to level 1)




5. Short press the on/off button to adjust the energy level 1-9.

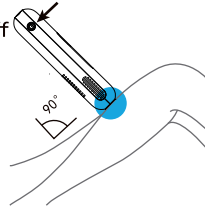


6. Hold the device vertically close to the skin so that it is pressed against the skin. After skin sensing the display on the screen changes to a constant flashing hair count, indicating that the epilator is ready for power and can be used

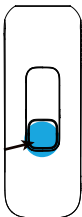


7. Press the flash button to emit the flash. Press and hold for 3 seconds to switch between auto mode and manual mode. Double-press the Flash button switching into HR/RA/SC mode.

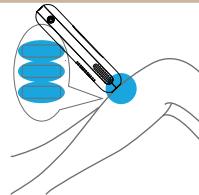

8. When using the hair removal device slide to the next area that needs to be irradiated onto the skin. Increase the level step by step within the recommended range and do not flash the same skin area multiple times.

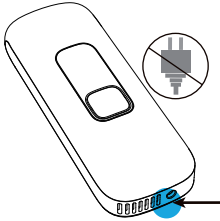

9. Long press the "power" button for 2 seconds to turn off the device. (Never directly plug off the plug when not power off, avoid it caused damages)

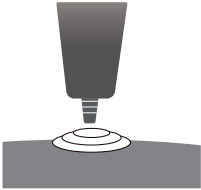

10. Apply skincare product to protect the skin. (Note: Avoid exposure to the sun within 24 hours after hair removal.)

## Applicable skin color and hair color


Select the suitable intensity level  
Please refer to the table to select the intensity level that suits your skin tone and hair color and confirm whether the product is suitable for you.

Skin Tone Hair Color	●	●	●	●	●	●
White/Grey	✗	✗	✗	✗	✗	✗
Light Blonde	✗	✗	✗	✗	✗	✗
Red	✗	✗	✗	✗	✗	✗
Dark Blonde Light Brown	✓	✓	✓	✓	✗	✗
Brown	✓	✓	✓	✓	✗	✗
Dark Brown	✓	✓	✓	✓	✗	✗
Black	✓	✓	✓	✓	✗	✗

## Applicable Parts and Recommended Use Time

**FOR WOMEN:** Head to toe usage including upper lip, chin, underarms, legs, arms, and bikini area.

**FOR MEN:** the ideal tool to tackle chest, back, arms, stomach, and leg hair.

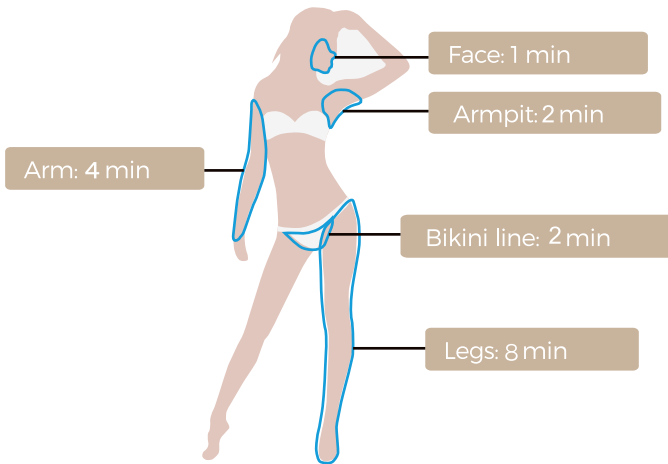

**Note:** The above recommended hair removal times are the result of research and testing. Individual hair removal times may vary from person to person.

## Common Problems and Solutions

Problems	Reason	Solution
The adapter or hair removal device heats up during use.	Normal	No action required.
The device is turned on, but the intensity level cannot be increased or decreased.	The device needs to be reset.	To restart the hair removal device, simply turn off the power and then turn it back on. Please contact the seller if you still can't adjust the light intensity.



Problems	Reason	Solution
After the instrument is powered on, it can't emit light normally when it is placed on the skin.	The instrument don't fit the skin vertically at 90°	After the instrument is vertically attached to the skin at 90°, the numbers on the LED screen flicker, indicating that the instrument is ready to work and can be used with normal flash.
After a long time of use, the numbers on the LED screen suddenly don't flicker, but the cooling fan is running normally.	The overheating protection system has been activated, suspending the flash work, while the fan is running normally to dissipate heat.	Don't turn off the device, let it cool down for a minute before continuing to use it. If it still doesn't flash after cooling down, try restarting.
The numbers on the LED screen flicker, but no flash is emitted when the flash button is pressed	Equipment system failure	reboot device. If it still does not flash after the second turn on, please contact the seller for after-sales service.
Pain feeling intolerable when using this product.	You have not remove the hair before using.	If it is still has hair at the position where need remove hair then these hair maybe burn when using this product then generate abnormal smell. Please remove the hair where need remove hair before using.
	The selected light energy too high.	The first time you use it is recommended to start from 1 gear, gradually strengthen to a bearable gear.
	The position where Can't use this device.	Can't use this product on Private parts, such as vagina, nipple or scrotum.
	This device unsuitable for your hair color or skin color.	IPL hair removal device is ideal for lighter skin tones and brown-to-dark hair. NOT applicable to black skin. NOT applicable to blonde, red, gray, white hair.
The position change to be red after hair removal.	Tiny red is normal appearance, harmless to human body, and will degrade soon .	No action required.

Problems	Reason	Solution
No result after using	The selected light strength much lower.	Execute the skin test to confirm whether you can use more higher light strength.
	You will be more likely to get results based on the following experimental suggestions:	1-4 weeks: 2 times every weeks; 5-8 weeks: 1 time every weeks; After 8 weeks: 1-2 times every months.
	No shaving before using	Shave off the hair on the surface of the skin every treatment to get better result.
The product unable to work	This product exist faults	Please contact the seller to handle it. The seller promises a 12-months unconditional warranty for any quality problems or no result issues no matter the Amazon return window is closed or not.

## QA

Question	Answer
30 days Amazon return window has passed. What can I do if my laser hair removal order has expired the return date?	Indeed, the 30 days isn't quite enough to effectively try the IPL hair removal. The seller hosting 12-months unconditional warranty. You could send message to the seller Via: ( Find "Your Orders" - find your this order ID - Click "Contact Seller").
Is this product suitable for pregnant women and children?	Since the hair removal device has not been tested in pregnant women and minors. It is not recommended to use by pregnant women and children.

Question	Answer
Does this laser hair removal device work on blonde hair ?	IPL hair removal device sends light pulses that are optimized to target the melanin in the hair. The light is absorbed by the melanin, which heats up, preventing unwanted hair from reaching the skin surface with continued use. Without melanin, this is impossible. Thus IPL does not work for blonde, gray and white hairs.
Do I need to shave the unwanted hair every treatment, or only for the first time?	This IPL device come with a razor. In order to get better result, you can shave the unwanted hair every treatment with this laser hair removal, not only the first time.
Is it ok to use over tattoos?	IPL hair removal can not use over tattoos.
How often should I have treatment?	You can follow the advice as follow: 1-4 weeks: 2 times every weeks; 5-8 weeks: 1 time every weeks; After 8 weeks: 1-2 times every months.
Is this safe to treat the entire bikini area with this? Full Brazilian	You can use laser hair removal permanent for bikini line. We don't recommend treating the genital area but consider completing a full bikini which will be very close, but not a full Brazilian.

